



Innovation in cryopreservation

◀ upstream ▶ downstream 🧩 single-use 📦 services

Advancing post-thaw cell recovery and growth with HyCryo and HyCryo-STEM cryopreservation media

Imagination at work.

HyCryo medium



Superior recovery for standard cell lines

Animal origin-free and serum-free cryopreservation medium for general use

With GE Healthcare HyCryo cryopreservation medium, you can store your samples with even more confidence, knowing that our validated media design advances the post-thaw recovery, viability, and growth of your cells.

HyCryo medium is a serum-free, animal origin-free product intended for cryogenic and preservation storage and biomanufacturing needs for your standard cell lines. With HyCryo medium, your precious cells are faithfully preserved

for your future research. Using HyCryo has been shown to increase recovery of viable cells versus results from other cryopreservation media alternatives.

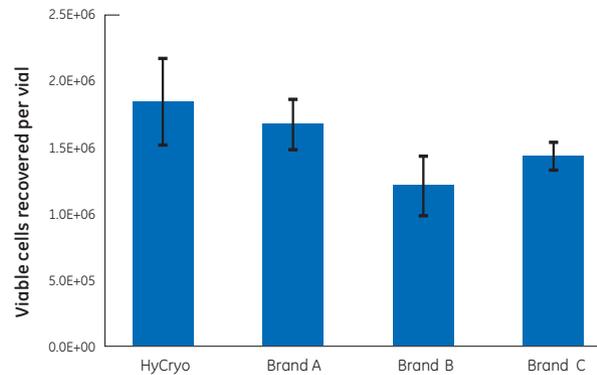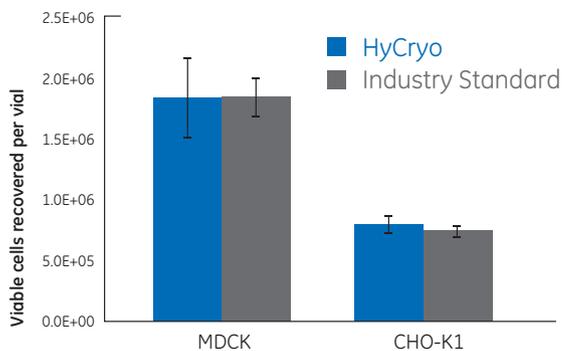

Ordering information

Product	Description	Size	Code no.
HyCryo	Cryopreservation medium for general use	100 mL	SR30001.02
HyCryo-STEM	Cryopreservation medium for stem cells	100 mL	SR30002.02




HyCryo-STEM medium

Designed specifically for stem cell storage


Serum-free formulation minimizes spontaneous differentiation

GE Healthcare HyCryo-STEM medium provides impressive post-thaw cell recovery and growth of pluripotent stem cells (ESC and iPSC) and neural progenitor cells such as mouse cortical stem cells (MCSCs) and rat cortical stem cells (RCSCs). The serum-free formulation maintains differentiation potential and minimizes spontaneous differentiation of stem cells.

Post-Thaw Recovery

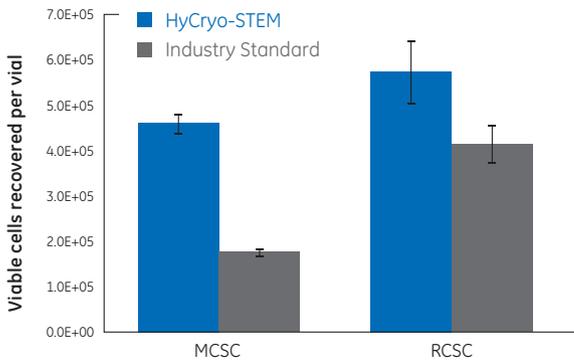


HyCryo-STEM medium was shown to enable greater recovery of viable neural progenitor stem cells compared to industry standard cryopreservation medium (growth medium containing 10% DMSO). Mouse cortical stem cells (MCSCs) were grown in adherent culture and rat cortical stem cells (RCSCs) were grown in neurosphere culture prior to freeze.

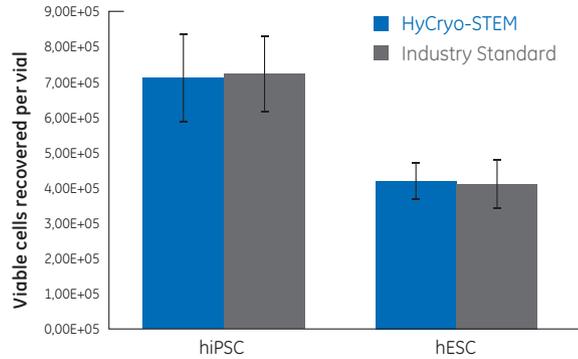


HyCryo-STEM medium achieves similar recovery of viable human stem cells compared to the industry standard cryopreservation medium (a cocktail of 25% serum replacement and 10% DMSO).

Post-thaw differentiation

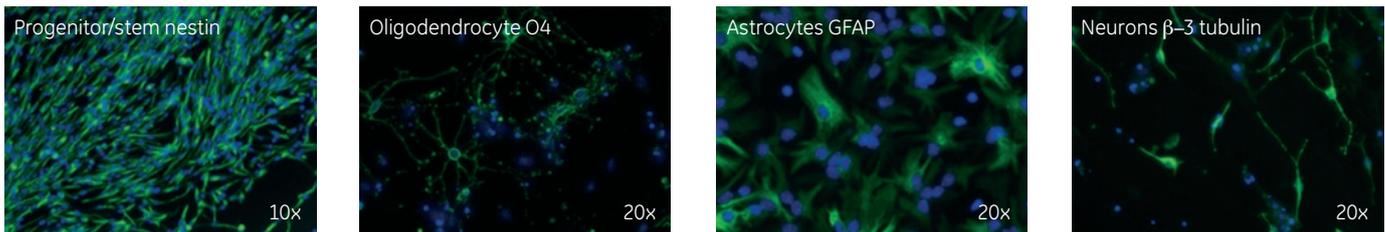

HyCryo-STEM medium maintains tripotency of neural progenitor stem cells post-thaw. Staining of mouse cortical stem cells grown as a monolayer and differentiated post-thaw demonstrate that cells frozen in HyCryo-STEM medium retain their stemness.

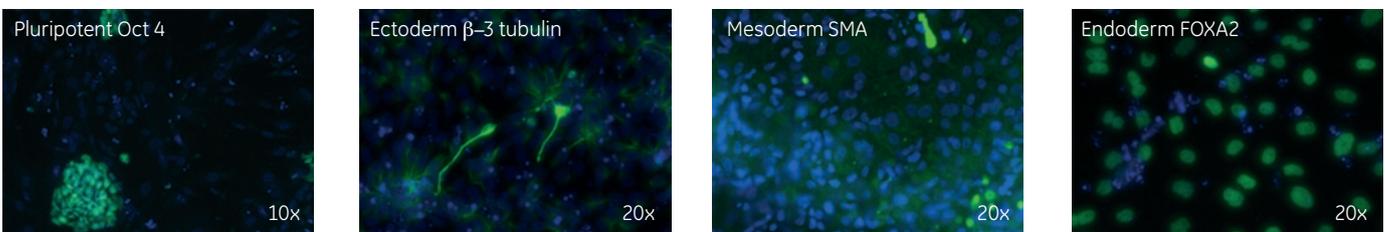

HyCryo-STEM medium maintains differentiation potential of human stem cells post-thaw. Human ES cells grown on MEFs post-thaw retained their pluripotency and were successfully differentiated into the three germ layers.

GE Healthcare UK Limited
Amersham Place, Little Chalfont
Buckinghamshire, HP7 9NA
UK

GE Healthcare Europe, GmbH
Munzinger Strasse 5, D-79111 Freiburg
Germany


GE Healthcare Bio-Sciences Corp.
800 Centennial Avenue
P.O. Box 1327
Piscataway, NJ 08855-1327
USA

GE Healthcare Japan Corporation
Sanken Bldg. 3-25-1, Hyakunincho
Shinjuku-ku, Tokyo 169-0073
Japan

For local office contact information, visit
www.gelifesciences.com/contact

GE Healthcare Bio-Sciences AB
Björkgatan 30
751 84 Uppsala
Sweden

www.gelifesciences.com/hyclone

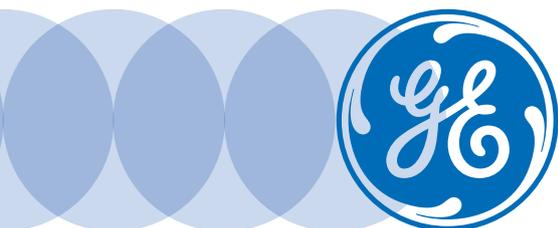

GE, imagination at work, and GE monogram are trademarks of General Electric Company.

HyClone is a trademark of General Electric Company or one of its subsidiaries.

© 2014 General Electric Company - All rights reserved.
First published Mar. 2014

All goods and services are sold subject to the terms and conditions of sale of the company within GE Healthcare which supplies them. A copy of these terms and conditions is available on request. Contact your local GE Healthcare representative for the most current information.